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
Ultrasonic Cleaner

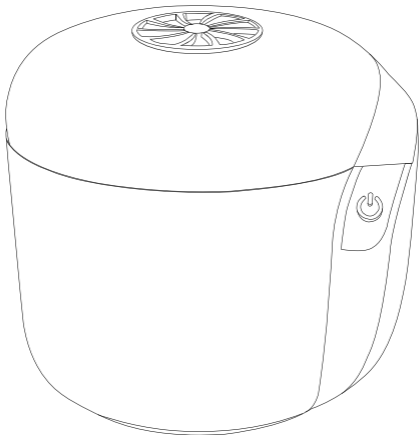
User Manual

MODEL: YBN-101

CONTENTS

English	- - - - -	01-03
Deutsch	- - - - -	04-06
Español	- - - - -	07-09
Français	- - - - -	10-12
Italiano	- - - - -	13-15

1. Product Overview

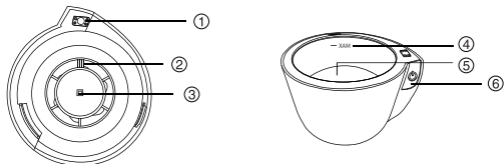1. Magnetic Lid
2. Fan
3. UVC Sterilization Lamp
4. Maximum Water Level Line
5. Cleaning Tank
6. Power On/Off


2. Specifications

Rated Input Voltage	12V (QC3.0)
Rated Current	2A
Rated Power	24W
Ultrasonic Cleaning Time	6 minutes (short press)
Air Drying Time	30 minutes (long press)
UVC Disinfection Time	6 minutes (Activated when both ultrasonic cleaning and air drying are in progress)
How to Shut Down	Short press the power button
Capacity	190ml


Note: Please use the provided power adapter to avoid insufficient power for cleaning or machine damage due to excessive power.

3. How to Use

a. Connect the power supply.

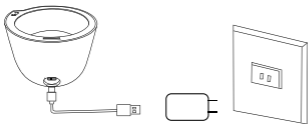b. Open the lid and place items in the cleaning tank.

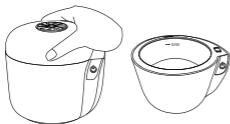c. Add water to the cleaning tank, ensuring it does not exceed the MAX water level line, preferably 1cm below.

d. Short press the power button for a 6-minute ultrasonic cleaning cycle; UV light operates with the closed lid during this time.

e. Long press the power button for a 30-minute air drying cycle; UV light runs for 6 minutes. The cleaner will auto-shutoff at the end of the cycle (Note: Empty water and wipe the tank clean before initiating the air drying cycle).

4. Item Cleaning Instructions

✓ **Items that can be cleaned:**

- **Oral care products:** retainers, braces, mouthguards, night guards, dentures, orthodontic appliances, aligners, clear aligner trays.
- **Jewelry:** gemstones, rings, toothbrush heads, earrings, necklaces.

× Items that cannot be cleaned:

- Items with adhesive gemstones: pearls, jade, emeralds, etc.
- Items with coatings or low-quality coatings
- Contact lenses

! Note:

- Tartar, plaque, stains, and tobacco stains in aligners cannot be completely cleaned.
- Discoloration and rust issues cannot be completely cleaned.

5. Daily Maintenance and Care

Clean with a soft sponge and rinse with water. For stubborn stains, add 1–2 drops of cleaning solution for better results; wipe dry with a cloth.



6. Accessories

- 1 x Charging Cable
- 1 x Power Adapter
- 1 x Ultrasonic Cleaner
- 1 x User Manual

7. Usage Precautions

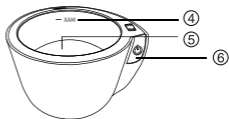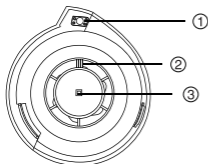- a. Keep out of reach of children to avoid burns.
- b. Add water before using the machine to prevent damage.
- c. Avoid direct eye contact with the UVC lamp to prevent discomfort.
- d. Do not fill water directly from the tap, or clean products under running water to prevent water from entering the power interface.
- e. Do not add disinfectants, sterilizers, or corrosive solvents to the water.

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