AGPTEK

Flossing & Toothbrush

User Manual


Model: ABB889

CONTENT

English	01–05
Deutsch	06–10
Français	11–15
Español	16–20
Italiano	21–25

Product Overview

1. Brush Head
2. Floss Nozzle
3. Storage Container
4. Water Tank
5. Hose
6. Leakproof Opening
7. Toothbrush On/Off /Mode Switch Button
8. Floss On/Off Button
9. Floss Speed Control Button
10. Flosser & Toothbrush Indicator Lights
11. Battery Power Indicator Lights
12. Holder
13. Type-C Charging Port
14. Charging Contacts

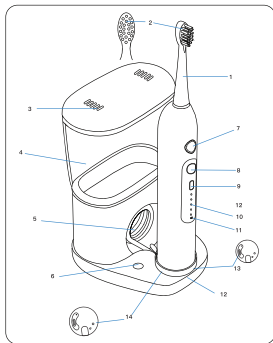Specifications

Material	Main Body: ABS+TPR; Hose: TPU.
	Water Tank: PC; Water Tank Holder: ABS.
Charging Port Input	5V/1A
Battery Capacity	1100mAh
Charging Time	<6H
Waterproof Rate	IPX7
Flossing Frequency	40–110PSI
Vibration Frequency	31200rpm
Indicator Light Description	Red Light: Brushing
	Blue Light: Flossing
	Purple Light: Flossing & Brushing

Package Contents


- 1 x Main Body
- 4 x Brush Heads
- 1 x Water Tank
- 1 x Hose
- 1 x Type-C Cable
- 1 x DC Cable
- 1 x Storage Bag
- 1 x User Manual

Instructions for Use

Steps to Use with Water Tank

1. Connect the brush head to the body

Place the brush head on the top of the body, and press it firmly until it snaps into place; to remove, pull the brush head hard to separate it from the body.

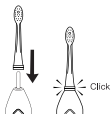
2. Connect the Hose to the Body

Place the body on the base and press the hose firmly into the water inlet on the back of the body.

3. Water Tank Installation


Take off the water tank cover, fill the water tank with water, and install the water tank back on the base and press it firmly.

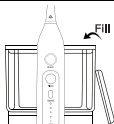
Steps to Use Without Water Tank

1. Connect the brush head to the body

Place the brush head on the top of the body, and press it firmly until it snaps into place; to remove, pull the brush head hard to separate it from the body.

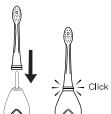2. Insert the hose into the body, insert the other side of the hose into a water container filled with water, and fix it with a fixing clip.

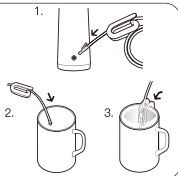
3 Modes

1. Floss

Aim the brush head between your teeth, bend over to the sink, close your lips slightly to prevent splashing, and press the FLOSS button.

Note: For first-time users and those with sensitive teeth, please use the slowest speed first, and then adjust the speed that suits you after you get used to it.

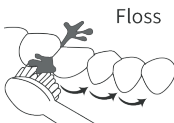2. Brush

Put the brush in your mouth, and press the BRUSH button.

3. Brush & Floss at the Same Time

Put the brush in your mouth, bend over to the sink, close your lips slightly, and press the BRUSH+FLOSS button.

Button Description

1. Brush

5 speeds are optional: Short press to power on and switch modes, and long press to power off.

2. Floss

Floss On/Off Button


3. Flossing Speed Adjustment Button

5 speeds are optional: From bottom to the top: 1 2 3 4 5.


30s Regional Reminder & 2mins Auto Power Off

The device pauses once every 30s to remind you to move the toothbrush to the next oral area. Clean the inside and outside of the teeth for 30s, repeat the operation on the other side, and clean all the teeth for a total of 2mins. The device will automatically power off after running for 2mins.

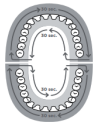Use of Mouthwash & Other Solutions

This device supports to floss with mouthwash and antibacterial solution. After using the solution, fill the water tank with warm water, power on the device and rinse the solution.

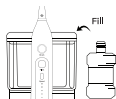Daily Care

1. After each use, turn off all functions, remove the toothbrush head, and rinse it with water.
2. Wipe the handle dry after each use.

Charging Methods

1. Type-C Charging.

2. Charging Contacts.

Battery Power Display

Red Light: Charging

Green Light: Fully Charged

Red Light Flashes: Low battery, 20% remain

Note: When using for the first time or not using for more than 6 months, the charging light may not light up for the first few minutes when charging starts.

Cleaning Care

Before cleaning, please disconnect the power.

General Cleaning/Removal of Deposits

1. Clean the product with a soft cloth and a mild, non-abrasive cleaner when necessary.
2. Deposits may build up in your device (depending on the mineral content of the water), which may affect product performance.
3. Clean the internal parts. Add 2 scoops of white vinegar to the warm water in the water tank, point the handle and brush head at the water tank; turn on the device and start working until the water in the water tank is used up, then rinse with clean warm water repeatedly.
4. Clean once every 1–3 months to ensure the best function.

FAQ

Problem	Reason	Solution
Nozzle cracking	The nozzle is badly damaged	Replace with a new brush head
Water tank leaking	Liquid remains in the tank	Keep the tank dry after use
	The water tank cover is not closed	Tighten the lid of the water tank
Water leakage at the joint between the main body and the nozzle	The nozzle is not connected properly	Remove the nozzle and reinstall
Liquid is still or under less pressure	The water tank is short of water	Fill the tank with water
	Working mode error	Switch to flossing mode

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